MEDFIRST.CARE


+33 5.32.63.00.14
medfirst.care
More Info !

AKED femoral stem extractor

Removing the femoral stem in revision hip surgery

In revision hip arthroplasty, this set is used once the prosthetic head has been removed, to free the femoral stem. A jaw extractor grips the implant taper, then a standard hammer fitted with an adjustable sliding weight performs the extraction.

The manufacturer intends it for stems with a neck carrying a modular head, whatever the model. Designed and made in France by AKED, it is distributed by Medfirst under reference 0081 01 00.

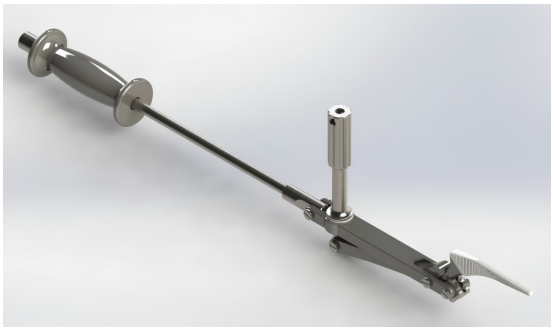

Technical specifications

Reference	0081 01 00
Contents	Femoral stem extractor and standard hammer with sliding weight
Materials	Stainless steel, polyethylene ring
Overall dimensions	65 × 12 × 6 cm
Use	Reusable, supplied non-sterile
UDI	03700751501944
Origin	France

Packaging and sterilization

The set is delivered non-sterile in an individual labelled plastic sleeve, with the instructions for use inside. As a reusable instrument, it is reprocessed according to the facility's cleaning and sterilization procedures, in line with these instructions.

CE marking

Class I medical device bearing the CE mark, covered by an EU declaration of conformity under Annex IV of Regulation (EU) 2017/745 (MDR). GMDN code 12696 and EMDN code L091101 (reusable orthopaedic prosthesis extractors).

Need help
or additional information ?

8 esplanade Compans Caffarelli - 31000 TOULOUSE, France
Email : ventes@medfirst.care

V1.0 - 10/11/2026

Set-up in four steps









1

Positioning

Once the head has been removed, the extractor jaws are placed around the stem taper.

2

Locking

The handle is tightened by hand until the taper is clamped between the jaws, which secures the instrument to the implant.

3

Axis adjustment

The hammer screws onto the extractor clevis, which is oriented in the chosen extraction direction.

4

Extraction

Blows delivered with the sliding weight gradually free the femoral stem.

Frequently asked questions

Which stems is it compatible with?

According to the manufacturer, all stems with a neck carrying a modular head.

What is a slide hammer?

A hammer whose weight slides along its shaft. Moving the weight generates the pulling impacts transmitted to the extractor, then to the stem.

Is it sterile or single use?

No. It is reusable and supplied non-sterile, to be reprocessed according to the manufacturer's instructions.

For other removal procedures, Medfirst also distributes a universal nail extraction kit and a broken screw extraction kit.

Need help
or additional information ?

8 esplanade Compans Caffarelli - 31000 TOULOUSE, France
Email : ventes@medfirst.care

V1.0 - 10/11/2026