


DRAIN'S KEHR - LATEX T-TUBE

PRESENTATION :

The Well Lead Latex T-tube is made of natural latex, coated with silicone. It is inserted into the common bile duct during surgery.

The Kehr's drain is used in open gallbladder surgery to drain bile. The rod is often connected to a bile bag.

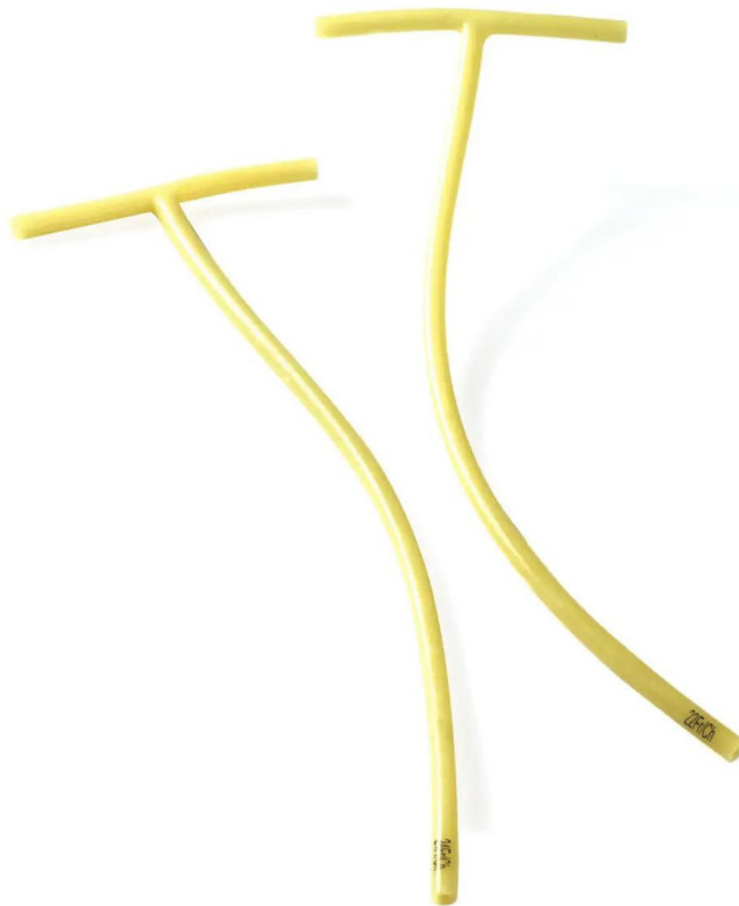

DESCRIPTION :

The Kehr's tube is a medical device used for the treatment of diseases of the gallbladder and bile ducts.

It is a specific type of drain used in surgery to drain bile from the gallbladder. It is often used during a cholecystectomy, a surgical procedure involving the removal of the biliary vesicle.

SIZES AVAILABLE :

Size	Length	Cat. No.	Size	Length	Cat. No.
10Fr	15cm x 30cm	N01D051020	10Fr	20cm x 40cm	N01D051030
12Fr	15cm x 30cm	N01D051220	12Fr	20cm x 40cm	N01D051230
14Fr	15cm x 30cm	N01D051420	14Fr	20cm x 40cm	N01D051430
16Fr	15cm x 30cm	N01D051620	16Fr	20cm x 40cm	N01D051630
18Fr	15cm x 30cm	N01D051820	18Fr	20cm x 40cm	N01D051830
20Fr	15cm x 30cm	N01D052020	20Fr	20cm x 40cm	N01D052030
22Fr	15cm x 30cm	N01D052220	22Fr	20cm x 40cm	N01D052230
24Fr	15cm x 30cm	N01D052420	24Fr	20cm x 40cm	N01D052430
26Fr	15cm x 30cm	N01D052620	26Fr	20cm x 40cm	N01D052630
28Fr	15cm x 30cm	N01D052820	28Fr	20cm x 40cm	N01D052830
30Fr	15cm x 30cm	N01D053020	30Fr	20cm x 40cm	N01D053030

CHARACTERISTICS :

It is shaped like a T, with one end positioned in the duct and the other end protruding through a small incision in the abdominal skin.

The drain, or T-tube, is inserted during the surgery, usually through the incision made in the abdominal wall.

Once inserted, the drain is guided to the gallbladder where it is secured in place with sutures.

The tube is then left in place for a variable period of time, ranging from a few days to several weeks, depending on the case.

While in place, the device allows for the drainage of bile that accumulates in the gallbladder after surgery.

This reduces pressure in the duct and helps to prevent complications such as leaks.

You can find our Kehr's drain, as well as additional information about this product, on our [dedicated page](#).

For more information, contact us

MEDFIRST

8 esplanade Compans Caffarelli 31000 TOULOUSE, France - RCS Toulouse B 831156641



+33 5.32.63.00.14



ventes@medfirst.care



www.medfirst.care